

Bench type, SD card real time data logger

LCR METER

Model : LCR-9185SD



Your purchase of this LCR METER with SD CARD DATA LOGGER marks a step forward for you into the field of precision measurement.

Although this meter a complex and delicate instrument, its durable structure will allow many years of use if proper operating techniques are developed. Please read the following instructions carefully and always keep this manual within easy reach.

OPERATION MANUAL

Caution Symbol




Caution :

- * Risk of electric shock !
- * During the measurement , do not open the cabinet.



Caution :

- * Do not apply the overload voltage, current to the input terminal !
- * Remove test leads before open the battery cover !
- * Cleaning - Only use the dry cloth to clean the plastic case !

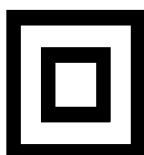

- * Double insulation

Environment Conditions

- * *Installation Categories III-1000V , IV-600V*
- * *Pollution Degree 2.*
- * *Altitude up to 2000 meters.*
- * *Indoor use.*
- * *Relative humidity 80% max.*

TABLE OF CONTENTS

1. FEATURES.....	1
2. SPECIFICATIONS.....	2
3. FRONT PANEL DESCRIPTION.....	9
4. MEASURING PROCEDURE.....	10
4-1. LCR Measurement	10
A). DCR Measurement.....	10
B). Z Measurement.....	10
C). C Measurement.....	11
a. Calibration.....	11
D).L Measurement.....	12
4-2 Data Hold.....	13
4-3 Data Record (Max., Min. reading).....	13
4-4 LCD Backlight ON/OFF.....	14
5. DATALOGGER.....	15
6. Download the saving data from the SD card to the computer (EXCEL software).....	19
7. ADVANCED SETTING	21
7-1 Set clock time (Year/Month/Date,Hour/Minute/ Second)..	21
7-2 Set sampling time (Seconds).....	22
7-3 Auto power OFF management.....	22
7-4 Set beeper sound ON/OFF.....	23
7-5 Decimal point of SD card setting.....	23
7-6 SD memory card Format.....	24
8. POWER SUPPLY from DC ADAPTER.....	25
9. BATTERY REPLACEMENT.....	25
10. FUSE REPLACEMENT.....	25
11. RS232 PC SERIAL OUTPUT.....	25
12. PATENT.....	27

1. FEATURES

- * 6000 counts ADC resolution.
- * High performance analog front end for impedance (Z) measurement.
- * Support Z / DCR measurement for LCR mode.
- * Four different test frequency are available :
100 Hz/120 Hz/1 KHz/10 KHz for L/C/R measurement.
- * Test AC signal level : 0.5 V rms typically.
- * Test range : (ex. F = 1 KHz)
L : 6.000 mH to 60.00 H
C : 6.000 nF to 600.0 uF
R : 600.0Ω to 20.00 MΩ
- * Min. source resistance : 120Ω typical.
- * 6 ratio resistor range used for LCR mode.
- * Real time Datalogger, save the into the SD memory card and can be downloaded to the Excel, extra software is no need.
- * Real time Datalogger, it Built-in Clock (year/month/date/hour/minute/second), sampling time set from 0 second to 3600 seconds.
- * Manual datalogger is available (set the sampling time to be 0 second).
- * Meet IEC 1010 CAT III 1000 V ,CAT IV 600 V safety requirement.
- * LCD with green light backligh, easy reading.
- * Accepts SD memory cards of up to 32 GB capacity.
- * Max. & Min. measurement value with recall.
- * Data hold.
- * Auto range with manual range selection.
- * Operates from 1.5 V (UM3/AA) x6 batteries or DC 9V Adapter.
- * Photo couple RS 232 computer serial interface.
- * Uses durable, long-lasting components, enclosed in strong, light weight ABS-plastic housing.

2. SPECIFICATIONS

2-1 General Specifications

Display	97 mm x 56 mm large LCD display	
Measurement	L/C/R Function selector frequency selector. SEL/PAL selector D/Q/θ selector	
Datalogger Sampling Time Setting range	Auto	0, 1, 2, 5, 10, 30, 60, 120, 300, 600, 1800, 3600 seconds
	Manual	Push the data logger button once will save data one time. @ Set the sampling time to 0 sec.
Data error no.	≤ 0.1% no. of total saved data typically.	
SD card Capacity	4 GB to 32 GB	
A/D counts no.	6000 counts.	
Range selection	Auto range with manual range selecting.	
Data hold	To freeze the display reading on the LCD display.	
Power On/Off management	Auto power off or manual power off. @ Details please refer page 21, 22	
Memory recall	Records Maximum & Minimum readings with recall.	
Data output	RS 232 PC serial interface, photo couple.	
Polarity	automatic switching, " - " indicates negative polarity.	
Zero adjustment	Automatic.	
Sampling time of display	Approx. 0.5 to 1 second.	
Operating Temp. & humidity	0 °C to 50 °C (32 °F to 122 °F), Max. 80% RH.	
Power supply	9 V (UM-3 1.5V X 6 PCS)	
	DC 9V adapter input (Linear 110V/220V) *AC/DC Power adapter is optional.	

Power consumption	DC 75mA Normal operation(w/o SD card save data and backlight off)
	When SD card save the data and LCD Backlight is OFF) will increase: Approx. DC 30 mA.
	<i>If LCD backlight on, the power consumption will increase approx.8 mA.</i>
Dimension	292 x 236 x 98 mm (11.5 x 9.3 x 3.9 inch).
Weight	1972 g/4.30 LB (w.o battery).
Accessories Included	Alligator clips(LN-TL03)..... 1 Set Instruction Manual..... 1 PC
Optional accessories	SMD test clip.....SMDC-21
	AC to DC 9V Adapter..... AA110V9 / AG220V9
	RS232 cable to USB connector..... USB-01
	SD memory card (4 G)..... SD-4GB
	Data Acquisition software..... SW-U801-WIN SW-E802

2-2 LCR General Specifications

Test frequency	100 Hz/120 Hz/1 KHz/10 KHz
MODE L/C/R	L/C/R Function selector Frequency selector D/Q/θ selector SER/PAL selector
Dissipation factor	0.000 to 9999
Quality factor	0.000 to 9999
θ measurement	± 90°
R measurement	0.00 Ω to 20.00 MΩ
Calibration	Open/Short calibration

2-4 LCR Electrical Specifications (23±5 °C)

Resistance (DCR)

Range	Accuracy	Remark
600 Ω	±(1.0% + 5d)	After calibration
6000Ω	±(1.0% + 5d)	
60 KΩ	±(1.0% + 5d)	
600 KΩ	±(1.0% + 5d)	
6000 kΩ	±(1.0% + 5d)	
20 MΩ	±(1.5% + 5d)	

Resistance(Z) (SER/PAL) 0.5V(rms)

Range	Accuracy	Accuracy	Remark
	100 Hz/120 Hz	1k Hz	
600 Ω	$\pm(1.2\% + 5d)$	$\pm(1.2\% + 5d)$	After calibration
6 K Ω	$\pm(1.2\% + 5d)$	$\pm(1.2\% + 5d)$	
60 K Ω	$\pm(1.2\% + 5d)$	$\pm(1.2\% + 5d)$	
600 K Ω	-----	$\pm(1.2\% + 5d)$	
6.000 M Ω	-----	$\pm(1.2\% + 5d)$	
20 M Ω	-----	$\pm(2.0\% + 5d)$	After calibration

Range	Accuracy	Remark
	10 kHz	
600 Ω	$\pm(1.2\% + 5d)$	After calibration
6 K Ω	$\pm(1.2\% + 5d)$	
60 K Ω	$\pm(1.2\% + 5d)$	
600 K Ω	-----	
6.000 M Ω	-----	
20 M Ω	-----	

- Remark :**
- * *All specifications are under in battery operation.*
 - * *If the equivalent measurement value is unstable, it is recommended to connect the GUARD terminal to the earth to obtain a stable reading.*
 - * *Don't apply voltage larger than 30 V to input terminals.*
 - * *When Auto range Function, test frequency default 1 KHz, Frequency button is disable.*

Capacitance (SER/PAL) : $D \leq 0.1$, 0.5V(rms)

Range	Accuracy	Accuracy	Remark
	100 Hz	120 Hz	
600 pF	-----	-----	
6 nF	-----	-----	
60 nF	$\pm(2.0\% + 5d)$	$\pm(2.0\% + 5d)$	
600 nF	$\pm(2.0\% + 5d)$	$\pm(2.0\% + 5d)$	
6 μ F	$\pm(1.5\% + 5d)$	$\pm(1.5\% + 5d)$	
60 μ F	$\pm(1.5\% + 5d)$	$\pm(1.5\% + 5d)$	
600 μ F	$\pm(1.5\% + 5d)$	$\pm(1.5\% + 5d)$	
6000 μ F	$\pm(2.5\% + 5d)$	$\pm(2.5\% + 5d)$	After calibration

Range	Accuracy	Accuracy	Remark
	1k Hz	10 kHz	
600 pF	-----	$\pm(2.0\% + 5d)$	After calibration
6 nF	$\pm(2.0\% + 5d)$	$\pm(1.5\% + 5d)$	After calibration
60 nF	$\pm(2.0\% + 5d)$	$\pm(1.5\% + 5d)$	
600 nF	$\pm(1.5\% + 5d)$	$\pm(1.5\% + 5d)$	
6 μ F	$\pm(1.5\% + 5d)$	$\pm(1.5\% + 5d)$	
60 μ F	$\pm(1.5\% + 5d)$	$\pm(2.5\% + 5d)$	
600 μ F	$\pm(2.5\% + 5d)$	-----	
6000 μ F	-----	-----	

Remark :

- * **All specifications are under in battery operation.**
- * **If the equivalent measurement value is unstable, it is recommended to connect the GUARD terminal to the earth to obtain a stable reading.**
- * **Don't apply voltage larger than 30 V to input terminals.**
- * **Discharge capacitor before measurement.**
- * **Precise measurement of surface mount components through optional test pieces.**
(SMD test clip, SMDC-21.)

Inductance (SER/PAL) : $D \leq 0.1$, 0.5V(rms)

Range	Accuracy	Accuracy	Remark
	100 Hz	120 Hz	
600 uH	-----	-----	
6 mH	-----	-----	
60 mH	$\pm(2.0\% + 5d)$	$\pm(2.0\% + 5d)$	
600 mH	$\pm(1.5\% + 5d)$	$\pm(1.5\% + 5d)$	
6 H	$\pm(1.5\% + 5d)$	$\pm(1.5\% + 5d)$	
60 H	$\pm(1.5\% + 5d)$	$\pm(1.5\% + 5d)$	After calibration
200 H	$\pm(2.5\% + 5d)$	$\pm(2.5\% + 5d)$	After calibration
Range	Accuracy	Accuracy	Remark
	1k Hz	10 kHz	
600 uH	-----	$\pm(2.5\% + 5d)$	After calibration
6 mH	$\pm(2.0\% + 5d)$	$\pm(2.0\% + 5d)$	
60 mH	$\pm(1.5\% + 5d)$	$\pm(1.5\% + 5d)$	
600 mH	$\pm(1.5\% + 5d)$	$\pm(1.5\% + 5d)$	
6 H	$\pm(1.5\% + 5d)$	$\pm(1.5\% + 5d)$	
60 H	$\pm(2.5\% + 5d)$	-----	After calibration
200 H	-----	-----	

Remark :

- * ***All specifications are under in battery operation.***
- * ***Discharge inductor before measurement.***
- * ***If the equivalent measurement value is unstable, it is recommended to connect the GUARD terminal to the earth to obtain a stable reading.***
- * ***Don't apply voltage larger than 30 V to input terminals.***
- * ***Precise measurement of surface mount components through optional test pieces.***
(SMD test clip, SMDC-21.)

2-5 LCR Scale range configuration

LCR mode			
Function mode	Frequency	Meas. Range	Min. resolution
Inductance (SER/PAL)	100/120Hz	60.00mH~200.0H	0.01mH
	1kHz	6.000mH~60.00H	0.001mH
	10kHz	600.0uH~6.000H	0.1uH
Capacitance (SER/PAL)	100/120Hz	60.00nF~6000uF	0.01nF
	1kHz	6.000nF~600.0uF	0.001nF
	10kHz	600.0pF~60.00uF	0.1pF
Resistance (SER/PAL)	100/120Hz	600.0 Ω ~60.00K Ω	0.1 Ω
	1kHz	600.0 Ω ~20.00M Ω	0.1 Ω
	10kHz	600.0 Ω ~60.00K Ω	0.1 Ω

3. FRONT PANEL DESCRIPTION

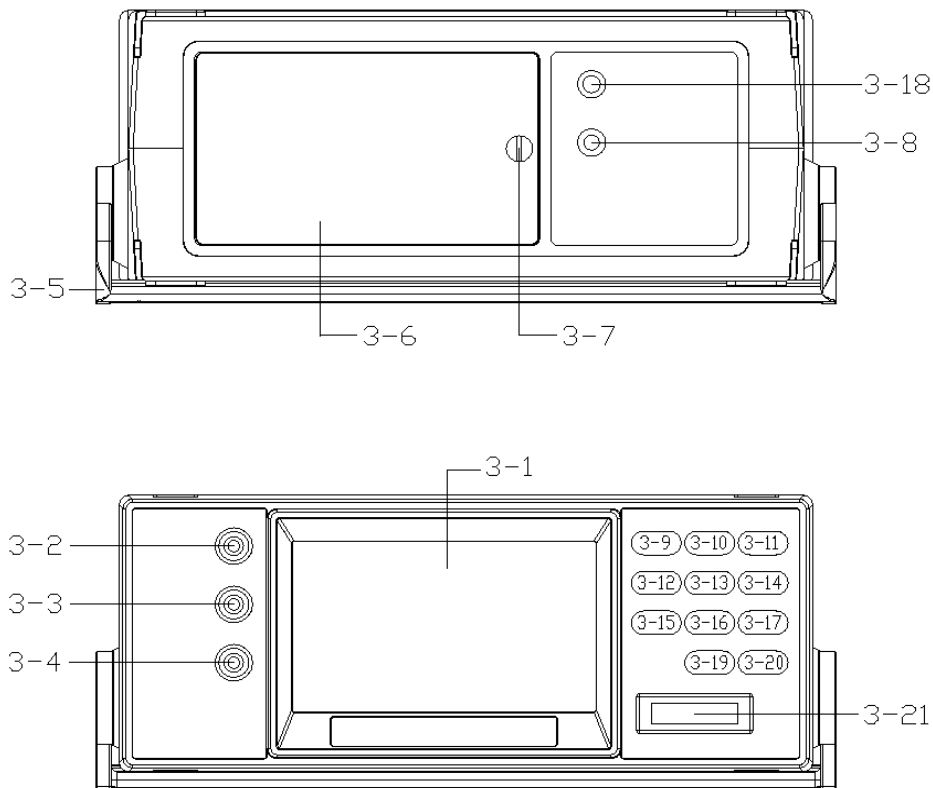


Fig. 1

- | | |
|-------------------------------|-----------------------------|
| 3-1 LCD Display | 3-12 FUNCTION Button |
| 3-2 LCR+ Input Terminal | 3-13 REC/▲ Button |
| 3-3 LCR- Input Terminal | 3-14 V/A/Hz Button |
| 3-4 GUARD Terminal | 3-15 RANGE Button |
| 3-5 Stand | 3-16 POWER/Backlight Button |
| 3-6 Battery Compartment/Cover | 3-17 TIME/SET Button |
| 3-7 Battery Cover Screws | 3-18 RS-232 Output Terminal |
| 3-8 DC 9V Power Adapter Input | 3-19 HOLD/▼ Button |
| 3-9 FREQUENCY Button | 3-20 ENTER/Log Button |
| 3-10 SER/PAR Button | 3-21 SD card Socket |
| 3-11 D/Q/θ Button | |

4. MEASUREMENT PROCEDURE

4-1 LCR Measurement

A). DCR Measurement

- 1) : Connect BLACK test lead into " - " terminal (3-3, Fig. 1).
- 2) : Connect RED test lead into " + " terminal (3-2, Fig. 1).
- 3) : Press the " Function Button " (3-12, Fig. 1) Select to the " DCR " Function.
- 4) : When LCD show the " AUTO " marker, the meter is under the " auto range " mode. Meter will select the suitable measurement range automatically.
- 5) : Under the operation of " auto range " mode, short press the " Range Button " (3-15, Fig.1) select to the Manual Range ,Manual Range can select 600Ω, 6KΩ, 60KΩ, 600KΩ, 6MΩ, 20MΩ, LCD will show the " MANU " marker, and Under the manual range operation, press and hold the " Range Button " (3-15, Fig.1) > 2 seconds, will return to auto range operation.
- 6) : When the connection is completed, the test value will be show on the LCD displayed.

B). Z Measurement

- 1) : Connect BLACK test lead into " - " terminal (3-3, Fig. 1).
- 2) : Connect RED test lead into " + " terminal (3-2, Fig. 1).
- 3) : Press the " Function Button " (3-12, Fig. 1) Select to the " Z " Function.
- 4) : When LCD show the " AUTO " marker, the meter is under the " auto range " mode. Meter will select the suitable measurement range automatically.
- 5) : Under the operation of " auto range " mode, short press the " Range Button " (3-15, Fig.1) select to the Manual Range ,Manual Range can select 600Ω, 6KΩ, 60KΩ, 600KΩ, 6MΩ, 20MΩ, LCD will show the " MANU " marker, and Under the manual range operation, press and hold the " Range Button " (3-15, Fig.1) > 2 seconds, will return to auto range operation.
- 6) : When "FREQ" key is pushed, the test frequency will be changed sequentially. There are four different test frequencies(100Hz/120Hz/1KHz/10KHz) can be selected. The LCR impedance scale range are depended on the test frequency, See the specifications.

- 7) : Use the " SER/PAL" key(3-10, Fig,1) short press to select the " serial " mode or "parallel" mode measurement.
- 8) : When the connection is completed, the test value will be show on the LCD displayed.

C). C Measurement

a. Calibration

In order to improve the accuracy for high / low impedance, it is recommended to do OPEN / SHORT calibration before the measurement.

OPEN / SHORT calibration procedures :

Press and hold the "▲ Button" (3-13, Fig, 1) and "▼ Button" (3-19, Fig, 1) over 2 seconds to start the open / short calibration , At the same time will be show text" CAL "and" OPEn " on the LCD panels , then press "ENTER Button"(3-20, Fig, 1) once , The LCD will show "ENT" and "30", press the "ENTER button" (3-20.fig, 1) once, At the same time 30 seconds will be countdown to zero. LCD will be show text "CAL" and "Short" on the LCD panels ,then shorting the LCR meter test input, and press "ENTER Button"(3-20, Fig, 1) once , The LCD will show "ENT" and "30", press the " ENTER button " once, At the same time 30 seconds will be countdown to zero. then the calibration data will be saved to EEPROM.

- 1) : Connect BLACK test lead into " - " terminal (3-3, Fig. 1).
- 2) : Connect RED test lead into " + " terminal (3-2, Fig. 1).
- 3) : Press the " Function Button " (3-12, Fig. 1) Select to the " C" Function.
- 4) : When LCD show the " AUTO " marker, the meter is under the " auto range " mode. Meter will select the suitable measurement range automatically.

- 5) : Under the operation of " auto range " mode, short press the " Range Button " select to the Manual Range ,Manual Range can select 600p,6nF,60nF,600nF, 6uF,60uF,600uF,6000uF .LCD will show the " MANU " marker, and Under the manual range operation, press and hold the " Range Button " (3-15, Fig.1) > 2 seconds, will return to auto range operation.
- 6) : When "FREQ" key is pushed, the test frequency will be changed sequentially. There are four different test frequencies(100Hz/120Hz/ 1KHz/10KHz) can be selected. The LCR impedance scale range are depended on the test frequency, See the specifications.
- 7) : Use the " D/Q/θ button " (3-11, Fig. 1) short press to select the " D " or "Q "measurement, or "θ "measurement, or " R "measurement.
- 8) : Use the " SER/PAL" key(3-10, Fig,1) short press to select the " serial " mode or "parallel" mode measurement.
- 9) : When the connection is completed, the test value will be show on the LCD displayed.

D).L Measurement

- 1) : Connect BLACK test lead into " - " terminal (3-3, Fig. 1).
- 2) : Connect RED test lead into " + " terminal (3-2, Fig. 1).
- 3) : Press the " Function Button " (3-12, Fig. 1) Select to the " L " Function.
- 4) : When LCD show the " AUTO " marker, the meter is under the " auto range " mode. Meter will select the suitable measurement range automatically.
- 5) : Under the operation of " auto range " mode, short press the " Range Button " select to the Manual Range ,Manual Range can select 6mH,60mH,600mH, 6H,60H,200H .LCD will show the " MANU " marker, and Under the manual range operation, press and hold the " Range Button " (3-15, Fig.1) > 2 seconds, will return to auto range operation.

Remark :

※ 200H range is limited to 100Hz and 120Hz measurement mode.

- 6) : When "FREQ" key is pushed, the test frequency will be changed sequentially. There are four different test frequencies(100Hz/120Hz/1KHz/10KHz) can be selected. The LCR impedance scale range are depended on the test frequency, See the specifications.
- 7) : Use the " D/Q/θ " button " (3-11, Fig. 1) short press to select the " D " or "Q "measurement, or "θ "measurement, or " R "measurement.
- 8) : Use the " SER/PAL" key(3-10, Fig,1) short press to select the " serial " mode or "parallel" mode measurement.
- 9) : When the connection is completed, the test value will be show on the LCD displayed.

4-2 Data Hold

During the measurement, press the " Hold Button " (3-19, Fig. 1) once will hold the measured value & the LCD will display a " HOLD " symbol. Press the " Hold Button " once again will release the data hold function.

4-3 Data Record (Max., Min. reading)

- 1) The data record function records the maximum and minimum readings. Press the " REC Button " (3-13, Fig. 1) once to start the Data Record function and there will be a " REC " symbol on the display.
- 2) With the " REC " symbol on the display :
 - a) Press the " REC Button " (3-13, Fig. 1) once, the " REC MAX " symbol along with the maximum value will appear on the display.
 - b) Press the " REC Button " (3-13, Fig. 1) again, the " REC MIN " symbol along with the minimum value will appear on the display.

- c) Press the " REC Button " (3-13, Fig. 1) again, the " MIN " symbol will be disappeared, LCD only show the " REC " symbol, the meter start the Data Record function again.
- 3) To exit the memory record function, just press the "REC Button " (3-13, Fig. 1) > 3 seconds continuously. The display will revert to the current reading.

4-4 LCD Backlight ON/OFF

- 1) After power ON, the " LCD Backlight " will light automatically. During the measurement, Push the "Backlight Button " (3-16, Fig. 1) once , will turn OFF the " LCD Backlight ".
- 2) Push the " Backlight Button " once will turn ON the " LCD Backlight " again.

5. DATALOGGER

5-1 Preparation before execute datalogger function

a. Insert the SD card

Prepare a " SD memory card " (4G to 32G, optional), insert the SD card into the " SD card socket " (3-21, Fig. 1).


b. SD card Format

If SD card just the first time use into the meter, it recommend to make the " SD card Format " at first. please refer chapter 7(page 21) & chapter 7-6 (page 24).

c. Time setting

If the meter is used at first time, it should to adjust the clock time exactly, please refer chapter 7 & chapter 7-1 (page 21).

d. Decimal format setting

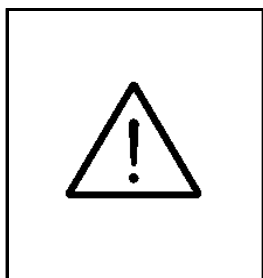

The numerical data structure of SD card is default used the " . " as the decimal, for example "20.6" "1000.53" . But in certain countries (Europe ...) is used the " , " as the decimal point, for example " 20, 6 " "1000,53". Under such situation, it should change the Decimal character at first, details of setting the Decimal point, refer to chapter 7 (page 21) & chapter 7-5 (page 23).

5-2 Auto Datalogger (Set sampling time ≥ 1 second)

a. Start the datalogger

Press the " LOG Button (3-20, Fig. 1) > 3 seconds continuously, the upper LCD will show the text of "Logger "then "Logger " symbol will flashing per sampling time, at the same time the measuring data along the time information will be saved into the memory circuit.

Remark :

- * *How to set the sampling time, refer to Chapter 7 (page 21) & Chapter 7-2 (page 22).*
- * *How to set the beeper sound is enable, refer to Chapter 7 (page 21) & Chapter 7-4 (page 23).*

b. Pause the datalogger

During execute the Datalogger function, if press the "LOG Button " (3-20, Fig. 1) once will pause the Datalogger function (stop to save the measuring data into the memory circuit temporally). In the same time the " Logger " symbol will stop flashing, the upper LCD will show the text of " Logger ".

Remark :

If press the " LOG Button " (3-20, Fig. 1) once again will execute the Datalogger again, the " Logger " symbol will flashing.

c.. Finish the Datalogger

During execute the Datalogger function, press the " LOG Button (3-20, Fig. 1) > 3 seconds continuously again will finish the Datalogger function, the " Logger " text will be disappeared and finish the Datalogger.

5-3 Manual Datalogger (Set sampling time = 0 second)

a. Set sampling time is to 0 second

Press the " LOG Button (3-20, Fig. 1) > 3 second , the LCD will show the " Position no. " , then press the " LOG Button " (3-20, Fig. 1) once, the " Logger " symbol will flashing once and Beeper will sound once, at the same time the measuring data along the time information will be saved into the memory circuit.

Remark :

During execute the Manual Datalogger, it can use the " ▲ Button " (3-13, Fig. 1) or " ▼ Button " (3-19, Fig. 1) to set the measuring position (1 to 99, for example room 1to room 99) to identify the measurement location , the Display will show P x (x = 1 to 99).

b. Finish the Datalogger

During execute the Datalogger function, press the " LOG Button (3-20, Fig. 1) > 3 seconds continuously again will finish the Datalogger function.

5-4 Check time information

During the normal measurement, If press " TIME Button " (3-17, Fig. 1) once , the lower LCD display will present the time information of Year/Month/Date , Hour/Minute/Second and the Sampling time information in sequence.

5-5 SD Card Data structure

- 1) When the first time, the SD card is used into the meter, the SD card will generate a folder :

LCA01

- 2) If the first time to execute the Datalogger, under the route LCA01\, will generate a new file name LCA01001.XLS. After exist the Datalogger, then execute again, the data will save to the LCA01001.XLS until data column reach to 30,000 columns, then will generate a new file, for example LCA01002.XLS
- 3) Under the folder LCA01\, if the total files more than 99 files, will generate a new route, such as LCA02\
- 4) The file's route structure :

LCA01\	LCA01001.XLS
	LCA01002.XLS

	LCA01099.XLS
LCA02\	LCA02001.XLS
	LCA02002.XLS

	LCA02099.XLS
LCAXX\

Remark :

XX : Max. value is 10.


6. Saving data from the SD card to the computer (EXCEL software)

- 1) After execute the Data Logger function, take away the SD card out from the " SD card socket " (3-21, Fig. 1).
- 2) Plug in the SD card into the Computer's SD card slot (if your computer build in this installation) or insert the SD card into the " SD card socket ". then connect the " SD card socket " into the computer.
- 3) Power ON the computer and run the " EXCEL software ".Down load the saving data file (for example the file name : LCA01001.XLS, LCA01002.XLS) from the SD card to the computer. The saving data will present into the EXCEL software screen (for example as following EXCEL data screens) , then user can use those EXCEL data to make the further Data or Graphic analysis usefully.

EXCEL data screen (for example)

Position	Date	Time	Ch1_Valu	Ch1_Unit	Ch2_Valu	Ch2_unit	Ch3_Valu	Ch3_unit
1	2017/6/6	15:03:49	47.3	nF	-89.7	Degree	1	KHz
2	2017/6/6	15:03:51	47.3	nF	-89.7	Degree	1	KHz
3	2017/6/6	15:03:53	47.3	nF	-89.7	Degree	1	KHz
4	2017/6/6	15:03:55	47.3	nF	-89.7	Degree	1	KHz
5	2017/6/6	15:03:57	47.3	nF	-89.7	Degree	1	KHz
6	2017/6/6	15:03:59	47.3	nF	-89	Degree	1	KHz
7	2017/6/6	15:04:01	57.72	nF	-89.7	Degree	1	KHz
8	2017/6/6	15:04:03	57.72	nF	-89.7	Degree	1	KHz
9	2017/6/6	15:04:05	57.72	nF	-89.7	Degree	1	KHz
10	2017/6/6	15:04:07	57.72	nF	-89.7	Degree	1	KHz
11	2017/6/6	15:04:09	57.73	nF	-89.7	Degree	1	KHz
12	2017/6/6	15:04:11	57.73	nF	-89.7	Degree	1	KHz
13	2017/6/6	15:04:13	68.2	nF	-89.7	Degree	1	KHz
14	2017/6/6	15:04:15	68.2	nF	-89.7	Degree	1	KHz
15	2017/6/6	15:04:17	68.2	nF	-89.7	Degree	1	KHz
16	2017/6/6	15:04:19	75.8	nF	-86.7	Degree	1	KHz
17	2017/6/6	15:04:21	78.7	nF	-89.7	Degree	1	KHz
18	2017/6/6	15:04:23	78.7	nF	-89.7	Degree	1	KHz
19	2017/6/6	15:04:25	78.7	nF	-89.7	Degree	1	KHz
20	2017/6/6	15:04:27	80.1	nF	-88.2	Degree	1	KHz
21	2017/6/6	15:04:29	89.3	nF	-89.7	Degree	1	KHz
22	2017/6/6	15:04:31	89.3	nF	-89.7	Degree	1	KHz
23	2017/6/6	15:04:33	89.3	nF	-89.7	Degree	1	KHz
24	2017/6/6	15:04:35	89.3	nF	-89.7	Degree	1	KHz
25	2017/6/6	15:04:37	99.8	nF	-89.7	Degree	1	KHz
26	2017/6/6	15:04:39	99.8	nF	-89.7	Degree	1	KHz
27	2017/6/6	15:04:41	99.8	nF	-89.7	Degree	1	KHz
28	2017/6/6	15:04:43	99.8	nF	-89.7	Degree	1	KHz

EXCEL graphic screen (for example)

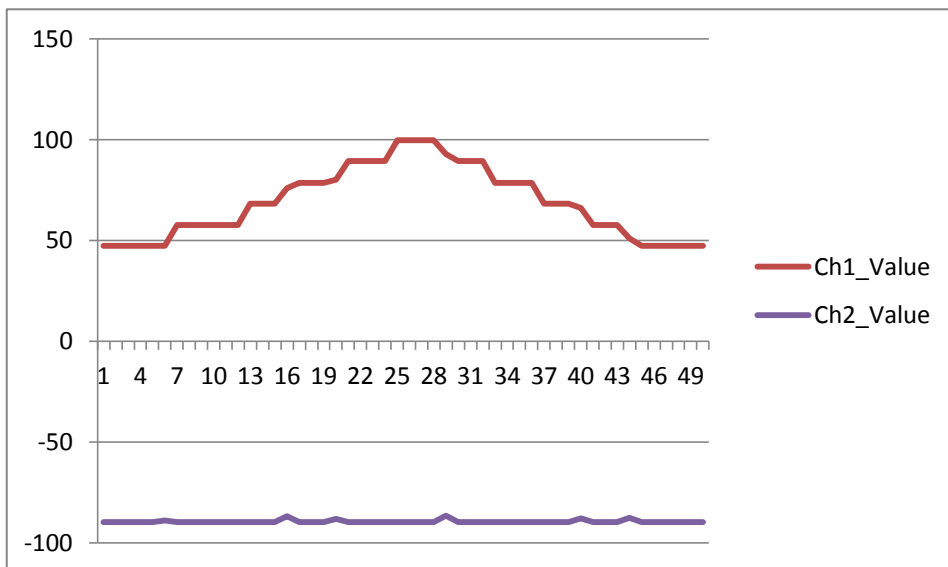

7. ADVANCED SETTING

Under do not execute the Datalogger function,
press the " SET Button " (3-17, Fig. 1) continuously at
least two seconds will enter the " Advanced Setting " mode.
then press the " SET Button " (3-17, Fig. 1) once a while
in sequence to select the six main function, the display
will show :

Sd F.....	SD memory card Format
dAtE.....	Set clock time (Year/Month/Date, Hour/Minute/ Second)
SP-t.....	Set sampling time
PoFF.....	Auto power OFF management
bEEP.....	Set beeper sound ON/OFF
dEC.....	Set SD card Decimal character

Remark :

*During execute the " Advanced Setting " function, if
all key button release > 10 seconds will exit the " Advanced Setting "
function, the LCD will return to normal screen.*

7-1 Set clock time (Year/Month/Date,Hour/Minute/ Second)

When the display show " dAtE "

- 1) Use the " ▲ Button " (3-13, Fig. 1) or " ▼ Button " (3-19, Fig. 1) to adjust the value (Setting start from Year value). After the desired value is set, press the " Enter Button " (3-20, Fig. 1) once will going to next value adjustment (for example, first setting value is Year then next to adjust Month, Date, Hour, Minute, Second value).

- 2) After set all the time value (Year, Month, Date, Hour, Minute, Second), press the " ENTER Button " (3-20, Fig. 1) once will save the time value.

Remark :

After the time value is setting, the internal clock will run precisely even Power off if the battery is under normal condition (No low battery power).

7-2 Set sampling time (Seconds)

When the display show " SP-t "

- 1) Use the " ▲ Button " (3-13, Fig. 1) or " ▼ Button " (3-19, Fig. 1) to adjust the value (0, 1, 2, 5, 10, 30, 60, 120, 300, 600, 1800,3600 seconds).
- 2) After the Sampling value is selected, press the " Enter Button " (3-20, Fig. 1) will save the setting function with default.

7-3 Auto power OFF management

When the lower display show " PoFF "

- 1) Use the " ▲ Button " (3-13, Fig. 1) or " ▼ Button " (3-19, Fig. 1) to select the upper value to " yES " or " no ".

yES - Auto Power Off management will enable.

no - Auto Power Off management will disable.

- 2) After select the upper text to " yES " or " no ", press the " ENTER Button " (3-20, Fig. 1) will save the setting function with default.

Remark :

Using the POWER ADAP Auto Power off function will be disable.

7-4 Set beeper sound ON/OFF

When the lower display show " bEEP "

- 1) Use the " ▲ Button " (3-13, Fig. 1) or " ▼ Button " (3-19, Fig. 1) to select the upper value to " yES " or " no ".

yES - Meter's beep sound will be ON with default.

no - Meter's beep sound will be OFF with default.

- 2) After select the upper text to " yES " or " no ", press the " ENTER Button " (3-20, Fig. 1) will save the setting function with default.

7-5 Decimal point of SD card setting

The numerical data structure of SD card is default used the " . " as the decimal, for example "20.6" "1000.53" . But in certain countries (Europe ...) is used the " , " as the decimal point, for example " 20,6 " "1000,53".

Under such situation, it should change the Decimal character at first.

When the lower display show " dEC "

- 1) Use the " ▲ Button " (3-13, Fig. 1) or " ▼ Button " (3-19, Fig. 1) to select the upper text to " USA " or " Euro ".

USA - Use " . " as the Decimal point with default.

Euro - Use " , " as the Decimal point with default.

- 2) After select the upper text to " USA " or " Euro ", press the " ENTER Button " (3-20, Fig. 1) will save the setting function with default.

7-6 SD memory card Format

When the lower display show " Sd F "

- 1) Use the " ▲ Button " (3-13, Fig. 1) or " ▼ Button " (3-19, Fig. 1) to select the upper value to " yES " or " no ".

yES - Intend to format the SD memory card


no - Not execute the SD memory card format

- 2) If select the upper to " yES ", press the " Enter Button " (3-20, Fig. 1) once again, the Display will show text " Sd F Ent " to confirm again, if make sure to do the SD memory card format, then press " Enter Button " once will format the SD memory clear all the existing data that already saving into the SD card.

8. POWER SUPPLY from DC ADAPTER

The meter also can supply the power supply from the DC 9V Power Adapter (optional). Insert the plug of Power Adapter into " DC 9V Power Adapter Input Socket " (3-8, Fig. 1).

9. BATTERY REPLACEMENT


- 1) When the left corner of LCD display show "  ", it is necessary to replace the battery. However, in-spec. measurement may still be made for several hours after low battery indicator appears before the instrument become inaccurate.
- 2) Loose the screw of the " Battery Cover Screws " (3-7, Fig. 1) and take away the " Battery Cover " (3-6, Fig. 1) from the instrument and remove the battery.
- 3) Replace with DC 1.5 V battery (UM3, AA, Alkaline/heavy duty) x 6 PCs, and reinstate the cover.
- 4) Make sure the battery cover is secured after changing the battery.

10. RS232 PC SERIAL INTERFACE

The instrument has RS232 PC serial interface via a 3.5mm terminal (3-18, Fig. 1).

The data output is a 16 digit stream which can be utilized for user's specific application.

A RS232 lead with the following connection will be required to link the instrument with the PC serial port.

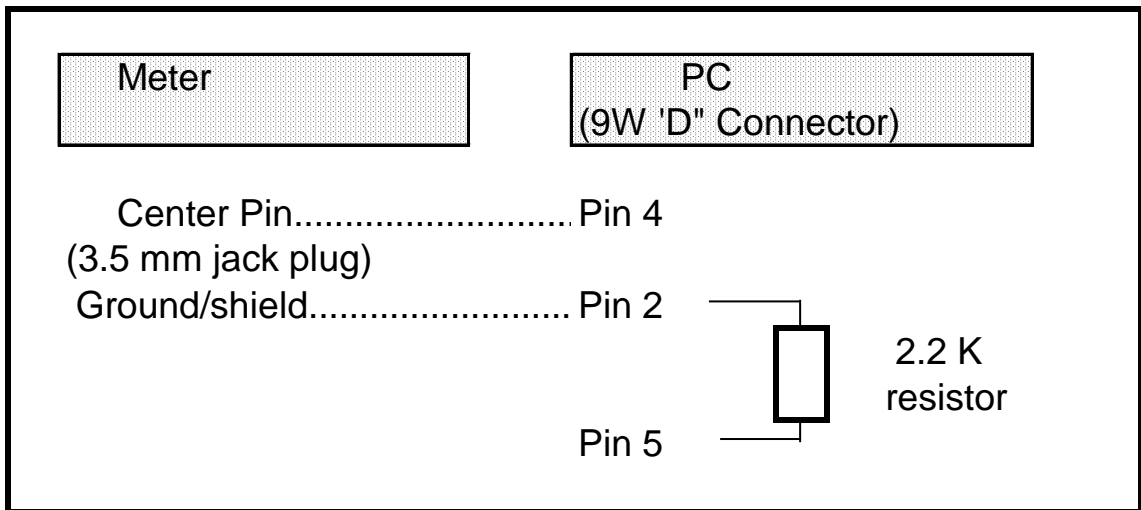

RS232 FORMAT : 9600, N, 8, 1

Baud rate	9600
Parity	No parity
Data bit no.	8 Data bits
Stop bit	1 Stop bit

The 16 digits data stream will be displayed in the following format :

D15 D14 D13 D12 D11 D10 D9 D8 D7 D6 D5 D4 D3 D2 D1 D0

Each digit indicate the following status :

D15	Start Word		
D14	4		
D13	1 ~ 2		
D12 & D11	Annunciator for Display		
	38 = OHM	39 = K OHM	40 = M OHM
	A4 = p F	43 = n F	44 = u F
	68 = u H	41 = m H	42 = H
	66 = m F		
D10	Polarity 0 = Positive 1 = Negative		
D9	Decimal Point(DP), position from right to the left 0 = No DP, 1= 1 DP, 2 = 2 DP, 3 = 3 DP		
D8 to D1	Display reading, D1 = LSD, D8 = MSD For example : If the display reading is 1234, then D8 to D1 is : 00001234		
D0	End Word		

12. PATENT

The meter (SD card structure) already get patent or patent pending in following countries :

Germany	Nr. 20 2008 016 337.4
JAPAN	3151214
TAIWAN	M 456490
CHINA	ZL 2008 2 0189918.5
	ZL 2008 2 0189917.0